



Ride® Designs

SEAT CUSHIONS &
BACK SOLUTIONS

Skin protection & position control



RIDE FORWARD

Improve postural control without compromise of skin protection

See clinical evidence on my-netti.com




FOR A DETAILED PRODUCT DEMONSTRATION, PLEASE CONTACT YOUR AREA MANAGER OR SEND US A MESSAGE TO info@my-netti.com

- Pre-shaped and anatomical design with great relief of the ischium and sacrum area
- Consistent pre-ischiatric stop to limit sliding
- Inner incontinence cover
- Breathable, non-compressable spacer fabric, to keep skin dry and cool
- Suitable for all ages, from adolescents to older people
- CAM wedges to correct asymmetrical pelvic challenges
- Polyurethane foam in the lower part and visco in the upper part
- Suitable for category 1 and 2 pressure injury
- For chairs with a seat width of 250-510 mm
- For chairs with a seat depth of 250-510 mm
- Other sizes available on request



More information:

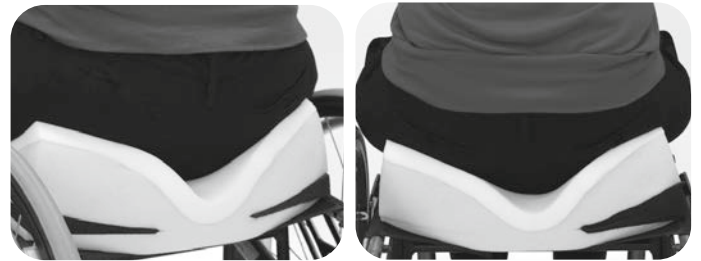


RIDE JAVA SEAT

Provides the highest level of skin protection and postural control

- Almost complete pressure relief of the high-risk areas due to the special design and combination of materials
- Suitable for users with moderate to a higher degree postural challenges because of the stabilising and pressure distributing features
- Weights about 30 % less than competitive cushions
- A vented base design, contoured top foam, and a breathable spacer mesh fabric cover work together to help keep skin cool and dry
- CAM wedges to correct asymmetrical pelvic challenges
- Medial and lateral thigh support available as option
- Suitable for category 3 and 4 pressure injury
- For chairs with a seat width of 360-510 mm
- For chairs with a seat depth of 360-510 mm
- Other sizes available on request



More information:

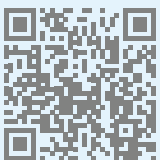

RIDE JAVA BACK

Adjustable Backrest with compact quick-release hardware

- Patented biomechanical design for optimal positioning of the pelvis and adapted seating comfort
- Biomechanical pelvic stop
- Compatible with Netti, MEYRA and TA wheelchairs
- Convenient adjustment of the back angle up to 20°
- Adjustable back adapter gives a 50 mm seat depth flexibility
- For extra support, fixed or swing-away lateral trunk support are available
- Available in seat width from 360-500 mm. Seat width can be adjusted by -25 mm or +50 mm
- Available in 3 heights: 350, 460 and 560 mm with 2 different depths of side supports (80 or 150 mm)

More information:



RIDE JAVA DECAF BACK

Pelvis and trunk support, comfort and enhanced mobility for children and small users

- Patented biomechanical design for child-friendly positioning
- Biomechanical pelvic stop
- Compatible with Netti, MEYRA and TA wheelchairs
- Back angle, height and depth adjustable
- Adjustable in transverse rotation to accompany trunk rotations
- For extra support, fixed or swing-away lateral trunk support are available
- Available in seat width 250-350 mm
- Different heights available from 250-380 mm depending on seat width
- Seat depth 50 or 76 mm

More information:





📍 MEYRA Group
Alu Rehabs ApS - part of MEYRA Group
Kløftehøj 8
DK-8680 Ry

✉ info@my-netti.com

Tel.: +45 878873 - 00

www.my-netti.com
